



TECH DATA SHEET

REPORTER VIRUS PARTICLES

DESCRIPTION

| | |
|--------------------------|--|
| Product | RVP-701L, SARS-CoV-2 Reporter Virus Particles (RVPs) |
| Lot | CL-112A |
| Strain | Wuhan-Hu-1 |
| Reporter | <i>Renilla luciferase</i> |
| Size | 1.0 mL/vial |
| Packaging | 20% FBS/DMEM |
| Recommended Input | 30uL per well (96-well plate) using the 293T-ACE2 stable cell line |
| Mycoplasma Test | Negative |
| Expiration Date | May 2024 |


SAFETY & HANDLING

| | |
|------------------------------|-------------------------------------|
| Shipping | Shipped on dry ice |
| Stability and Storage | Store at $\leq -80^{\circ}\text{C}$ |

SARS-CoV-2 RVPs are used to test the ability of serum, antibodies, and drugs to neutralize the infectivity of SARS-CoV-2 spike protein. RVPs display antigenically correct spike protein pseudotyped on replication-incompetent virus particles that contain a heterologous lentiviral (HIV) core. RVPs are capable of a single round of infection and carry a genome that expresses either a GFP or luciferase optical reporter gene upon infection. RVPs are produced in HEK-293T cells using three separate plasmids, encoding the spike protein, a lentiviral gag polyprotein, and a reporter gene. RVPs are created using a second-generation lentiviral system with components that are highly unlikely to recombine to produce a fully infectious virus (requiring 3 separate recombination events to do so). However, lentiviruses are capable of genomic integration and RVPs are derived from biological materials so should be handled with caution within a BSL2 or enhanced BSL2 laboratory environment. RVPs are not to be used in humans or in animals raised for food.

Thaw tubes in a 37°C water bath for 3 minutes and place on ice until ready to use. RVPs will appear as a translucent, pink solution. Gently mix prior to use and pulse tube for 3 seconds at high speed in a tabletop microfuge to recover all volume from the tube. Excessive vortexing of RVPs should be avoided. Re-freezing of RVPs is not recommended.

NEUTRALIZATION DATA

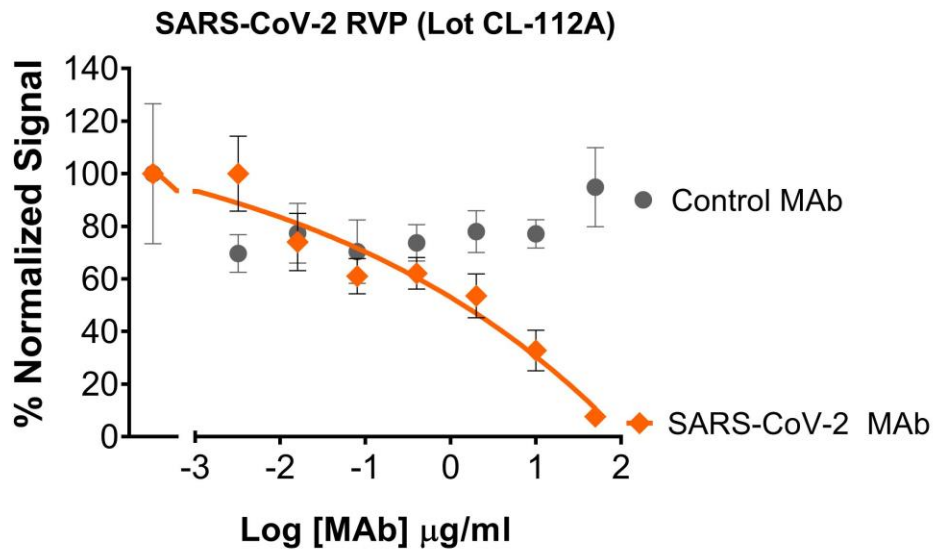

Infectivity determined in HEK-293T cells transiently over-expressing ACE2 and TMPRSS2 using 5 μl of SARS-CoV-2 RVPs in a 384-well plate. *Renilla* luciferase activity measured using the Promega *Renilla*-Glo luciferase assay system (Promega # E2710). Sample luminescence was read using a PerkinElmer Envision plate reader.