



TECH DATA SHEET

TiterSafe Hemagglutinating Particles

DESCRIPTION

Product	TiterSafe Hemagglutinating Particles, Influenza A
Catalog Number	HAP-1217
Lot	IAH-696A
Subtype	H5N1
Strain	dairy_cattle/Texas/24-008749-003-original/2024
Size	200 µL/vial
Packaging	DMEM
Dilution Ratio at 1 HA Unit	1:512
Dilution Ratio at 4 HA Units	1:128
Expiration Date	June 2026

SAFETY & HANDLING

Shipping	Shipped on dry ice
Stability and Storage	Store at ≤ -80°C upon receipt.


TiterSafe particles are optimized for use in HA and HAI assays. TiterSafe particles display antigenically correct HA/NA protein pseudotyped on replication-incompetent virus particles that contain a heterologous lentiviral (HIV) core. TiterSafe particles are capable of a single round of infection. TiterSafe particles are produced in HEK-293T cells using three separate plasmids, encoding the HA, NA, and lentiviral gag polyprotein.

TiterSafe particles are derived from biological materials, so they should be handled with caution within a laboratory environment. TiterSafe particles are not to be used in humans or animals raised for food.

Thaw tubes in a 37°C water bath for 1 minute and place on ice until ready to use. TiterSafe particles will appear as a translucent solution. Gently mix before use and pulse the tube for 3 seconds at high speed in a tabletop microfuge to recover all volume from the tube. Avoid vortexing TiterSafe particles. Re-freezing of TiterSafe particles is not recommended.

HA DATA

Lot IAH-696A HA assay

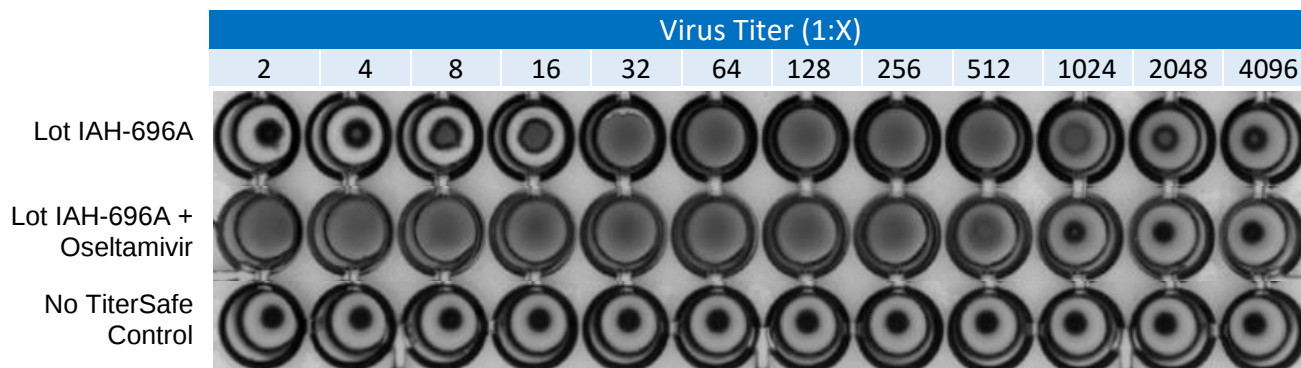

An HA assay was performed using 1% Guinea Pig RBCs in a 96-well plate. TiterSafe was treated with and without 1mg/ml Oseltamivir to counteract NA activity, which may interfere with the TiterSafe-RBC lattice at high TiterSafe concentrations. The HA assay was used to determine an Influenza A TiterSafe titer at 4 HA units.