



TECH DATA SHEET

REPORTER VIRUS PARTICLES

DESCRIPTION

Product	RVP-810G, Animal CoV Reporter Virus Particles (RVPs)
Lot	CG-861.4A
Species of Isolation	<i>Rhinolophus sinicus</i> (Horseshoe Bat)
Strain	WIV1
Reporter	GFP
Size	1.0 mL/vial
Packaging	20% FBS/DMEM
Viral Titer	5.30×10^6 TU/ml†
Recommended Input	5µL per well (96-well plate) for ≥ 20% infectivity in a flow assay*
Mycoplasma Test	Negative
Expiration Date	June 2028

SAFETY & HANDLING

Shipping	Shipped on dry ice
Stability and Storage	Store at ≤ -80°C upon receipt

* Determined in the 293T-ACE2 stable cell line

Animal CoV RVPs are used to test the ability of serum, antibodies, and drugs to neutralize the infectivity of Animal Coronaviruses. RVPs display antigenically correct spike protein pseudotyped on replication-incompetent virus particles that contain a heterologous lentiviral (HIV) core. RVPs are capable of a single round of infection and carry a genome that expresses either a GFP or luciferase optical reporter gene upon infection. RVPs are produced in HEK-293T cells using three separate plasmids, encoding the spike protein, a lentiviral gag polyprotein, and a reporter gene.

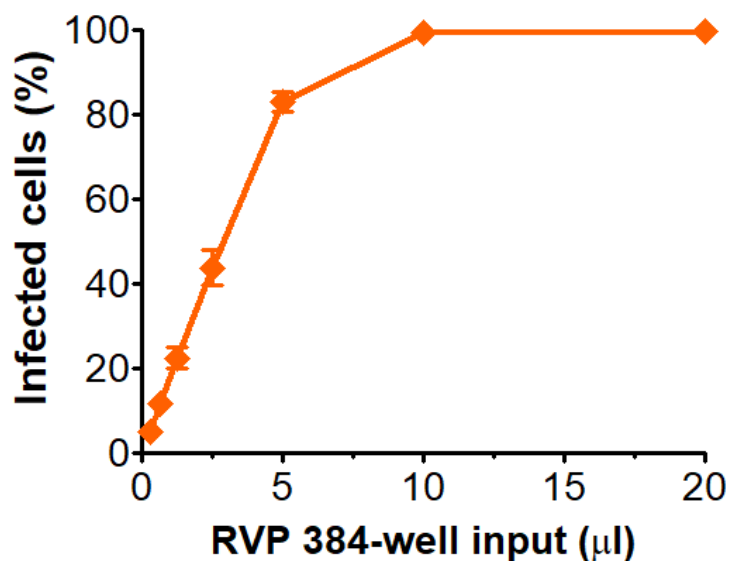
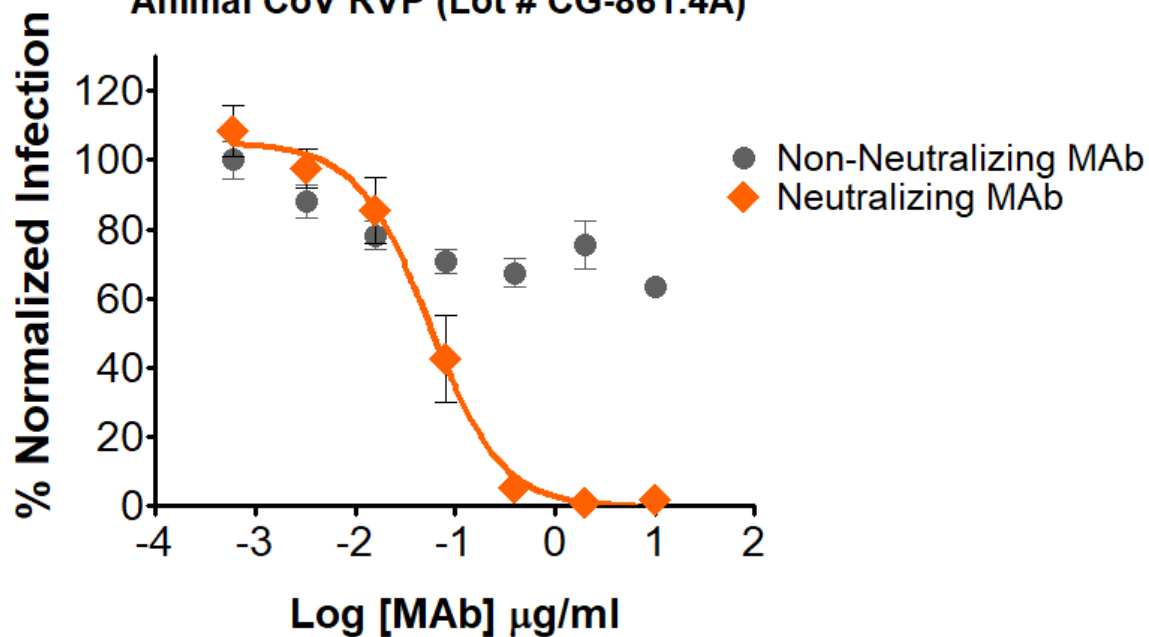
RVPs are created using a second-generation lentiviral system with components that are highly unlikely to recombine to produce a fully infectious virus (requiring 3 separate recombination events to do so). However, lentiviruses are capable of genomic integration and RVPs are derived from biological materials so should be handled with caution within a BSL2 or enhanced BSL2 laboratory environment. RVPs are not to be used in humans or in animals raised for food.

Thaw tubes in a 37°C water bath for 3 minutes and place on ice until ready to use. RVPs will appear as a translucent, pink solution. Gently mix prior to use and pulse tube for 3 seconds at high speed in a tabletop microfuge to recover all volume from the tube. Vortexing of RVPs should be avoided. Re-freezing of RVPs is not recommended.

† Titer is calculated using the estimated number of cells at the end of the experiment, using the following equation:

$$\left[\left(\frac{\% \text{ GFP positive cells}}{100} \right) \times \left(\frac{30,000 \text{ cells}}{1 \text{ well}} \right) \times \left(\frac{1 \text{ well}}{\text{RVP input } (\mu\text{L})} \right) \right] \times \frac{1000 \mu\text{L}}{1 \text{ mL}} = X \frac{\text{Transforming Units (TU)}}{\text{mL}}$$

INFECTIVITY AND NEUTRALIZATION DATA

Animal CoV RVPs (Lot # CG-861.4A)**Animal CoV RVP (Lot # CG-861.4A)**

Infectivity and neutralization determined in HEK-293T cells stably over-expressing ACE2 (Integral Cat# C-HA102). Infectivity data represents the average of three independent vials, each tested in quadruplicate.

Neutralization utilized 5 μl of Animal CoV RVPs in a 384-well plate. GFP positive cells were detected with an Intellicyt iQue flow cytometer using the BL-1 channel (Ex. 488 nm, Em. 530).